



**CPAPMAN.COM**  
BEST SERVICE SINCE 1997

*Click Here To Order Today*


# Sizing guide

## Instructions:

- 1 Print out this Adobe PDF File at 100% on your printer. Use the print accuracy check at right to ensure printing at the proper size.
- 2 Fold the bottom part of this page on the dashed line.
- 3 Cut along the dotted line around the fitting gauge.
- 4 Place the fitting gauge on your nose as shown in the picture below.
- 5 Determine the appropriate size by which area your nose reaches.  
(S – Small, M – Medium or L - Large)


Philips Respironics always recommends consulting a physician or registered Sleep Therapist in fitting a mask to ensure a proper sizing and fit.

## Print accuracy check

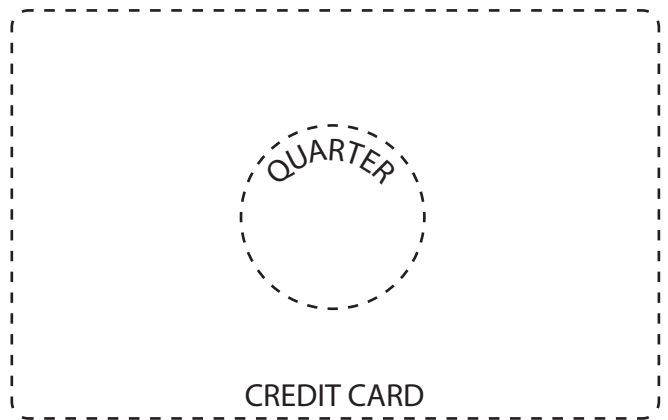

## Fitting gauge

CUT along dotted line

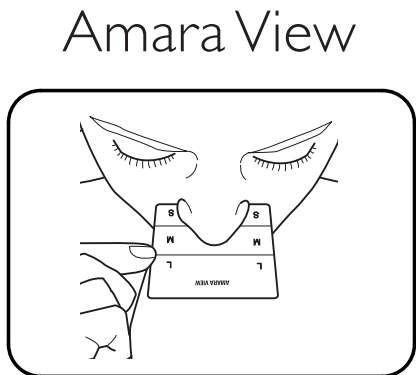

FOLD along dashed line

